



**savyonDIAGNOSTICS**

*member of the gamida diagnostics division*

**savvy•gen s**

## Genomic DNA Extraction Kit

**REF** 683-01; 96 extractions (96 wells plate)  
684-01; 48 extractions (48 wells plate)

Store at 15°C-35°C

For use with the Nextractor® NX-48S Instrument

For Professional Use Only **CE IVD**



**Savyon® Diagnostics Ltd.**  
3 Habosem St. Ashdod 7761003  
ISRAEL  
Tel.: +(972).8.8562920  
Fax: +(972).8.8523176  
E-mail: [support@savyondiagnostics.com](mailto:support@savyondiagnostics.com)



**European Authorized Representative:**  
**Obelis s.a.**  
Boulevard Général Wahis 53  
1030 Brussels, BELGIUM  
Tel: +(32) 2. 732.59.54  
Fax: +(32) 2.732.60.03  
E-Mail : [mail@obelis.net](mailto:mail@obelis.net)

## Table of Contents

|                                                  |          |
|--------------------------------------------------|----------|
| <b>Introduction .....</b>                        | <b>3</b> |
| Intended Use .....                               | 3        |
| Principle of Procedure .....                     | 3        |
| <b>Materials &amp; Equipment .....</b>           | <b>4</b> |
| Materials / Reagents Provided .....              | 4        |
| Additional Equipment and Material Required ..... | 4        |
| Kit Validity .....                               | 4        |
| <b>Warnings &amp; Precautions .....</b>          | <b>4</b> |
| <b>Experimental Process .....</b>                | <b>5</b> |
| Experimental Methods .....                       | 5        |
| Sample collection & Storage .....                | 5        |
| Test procedure.....                              | 5        |
| <b>Troubleshooting.....</b>                      | <b>8</b> |

# INTRODUCTION

## Intended Use

The Savvygen™ S Genomic DNA Extraction Kit is an automatic extraction system for the isolation of high quality DNA from Whole Blood, Buffy coat swab, Buccal swab, and saliva. This extraction kit is designed to use with Nextractor® NX-48S to provide high-yield and quality DNA from samples.

**For *in-vitro* professional diagnostic use.**

## Principle and procedure

The Savvygen™ Nucleic Acid Extraction technology is based on magnetic beads isolating of nucleic acid (NA) from Biological samples. The basic of this method is the use of Silica beads, capable of binding the NA in the presence of a chaotropic agent. This method is simple, rapid, and reliable method for the small-scale purification of NA from Biological sample.

The Savvygen™ Extraction kit is designed to use with Nextractor® system which offers fully automated extraction of up to 48 specimens at single and within 16 minutes.

The Savvygen™ Extraction kit procedure is performed in the following steps:

- 1) The specimens are placed into the sample well plate which contains in each well all the necessary reagents for the assay.
- 2) The well plate is placed on the loader tray of the Nextractor® system.
- 3) The user choose the desired protocol and start the program.

The extraction process for isolating nucleic acid from starting material by magnetic beads consist of the following 3 steps (Figure 1):

The first step of the extraction includes Lysis of the specimen in the presence of large amount of chaotropic substance and Binding of releasing-Nucleic Acid d to silica beads (magnetic bead).

The second step of extraction is washing silica beads. In this step, silica beads are washed in washing solution several times to remove contaminants such as PCR inhibitors.

The third step of extraction is DNA Elution. In this step, bonded nucleic acids is separated from the silica beads. Pure NA are eluted into buffer by decreasing the concentration of chaotropic substance

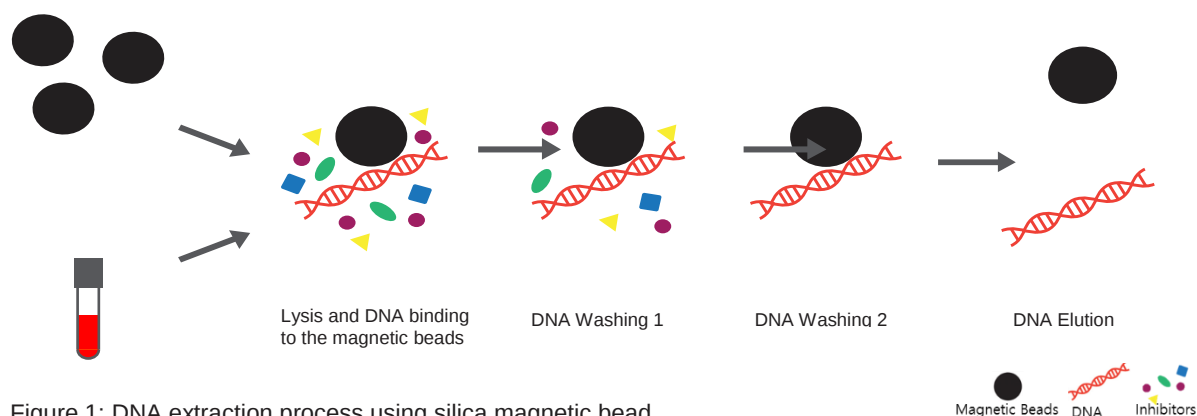


Figure 1: DNA extraction process using silica magnetic bead

# MATERIALS & EQUIPMENT

## Materials/ Reagents Provided

| Product Description                                                    | Contents                                             | Specimen type                                                 |
|------------------------------------------------------------------------|------------------------------------------------------|---------------------------------------------------------------|
| Savvygen™ S Genomic DNA Extraction Kit<br>Cat.# 683-01; 96 extractions | Extraction Plate (Genomic S)<br>(96 well x 4 plates) | Whole Blood,<br>Buffy coat swab,<br>Buccal swab and<br>saliva |
|                                                                        | Rods Plastic Strip x 12 units                        |                                                               |
|                                                                        | Proteinase K x 2 Tubes (2mL)                         |                                                               |
| Savvygen™ S Genomic DNA Extraction Kit<br>Cat.# 684-01; 48 extractions | Extraction Plate (Genomic S)<br>(32 well x 6 plates) |                                                               |
|                                                                        | Rods Plastic Strip x 6 units                         |                                                               |
|                                                                        | Proteinase K x 1 Tube (2mL)                          |                                                               |

## Additional Equipment and Material Required

- Disposable powder-free gloves
- Appropriate volume pipettes
- Sterilized, filtered pipette tips
- Vortex mixer
- Tabletop centrifuge
- Nextractor® (NX-48S, Genolution Inc. KOREA). Distributed by Savvyon Diagnostics Ltd. Please refer to the instrument manual for more details.
- Magnetic Bead Separation Rack, Magnetic Bead Holder

## Kit Validity

This kit is valid for 12 months when stored at 15-35°C.

## WARNINGS & PRECAUTIONS

- All samples must be treated as potential biohazards. Wear appropriate protective eyewear, clothing, and gloves.
- Avoid direct skin contact with kit reagents. In case of contact, wash immediately and thoroughly with water.
- Minimize the inhalation of chemicals. Do not leave chemical containers open.
- All work should be conducted in properly equipped facilities for safety reason (i.e. physical containment devices).
- Individuals should be trained according to the relevant regulation and requirements of the company/institutions prior to working with potentially infectious materials.

The reagent well plates contain Ethanol and Chaotropic Salt. These substances should be considered flammable, harmful and irritants. The Savvygen™ S Genomic DNA Extraction Kit and reagent well plates

are designed to be used with potentially infectious substances. Users should wear appropriate personal protective equipments (e.g. gloves and lab coat) when handling infectious substances.

#### **Extraction Plate (Genomic S)**

|                                                                                   |                                                                                                                                                                                                                                                      |
|-----------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  | <p><b>Harmful</b> – may cause sensitization skin contact. Avoid contact with skin, wear suitable gloves. Harmful to aquatic organizations may cause long-term adverse effects in the aquatic environments.<br/> <b>Contains:</b> Chaotropic Salt</p> |
|  | <p><b>Highly flammable</b>– keep away from any source of ignition, no smoking<br/> <b>Contains:</b> Ethanol</p>                                                                                                                                      |

## EXPERIMENTAL PROCESS

### **Experimental Methods**

- The Savvygen™ S Genomic DNA Extraction Kit should be used with Nextractor® NX-48S system
- All extractions should be processed in temperature between 15-35°C.
- Our extraction kit should be kept in the temperature between 15-35°C. Otherwise the result might be affected.

### **Sample collection & Storage**

#### Sample collection:

The Savvygen™ S Genomic DNA Extraction Kit is optimized for DNA extraction from various clinical samples as specified in the intended use. DNA should be isolated from clinical sample within one day after collection.

#### Sample storage:

The clinical samples are suggested to be stored for 1 day at 2-8°C. For longer storage, the clinical samples should be kept at -20°C. .

### **Test Procedure**

#### ***DNA extraction from Whole Blood, Buffy Coat swab, Buccal swab, saliva sample***

1. Determine the number of samples to be tested and carefully remove the sealing film of the Extraction plate.
2. Vortex samples for 10 sec.
3. Dispense 200µL of sample and 40µL of Proteinase K into the samples wells (in a 96 wells plate: apply samples to the 1<sup>st</sup>, 5<sup>th</sup>, and 9<sup>th</sup> columns; in a 32 wells plate: apply samples to the 1<sup>st</sup> column) of the figure extraction plate (See figure 2).

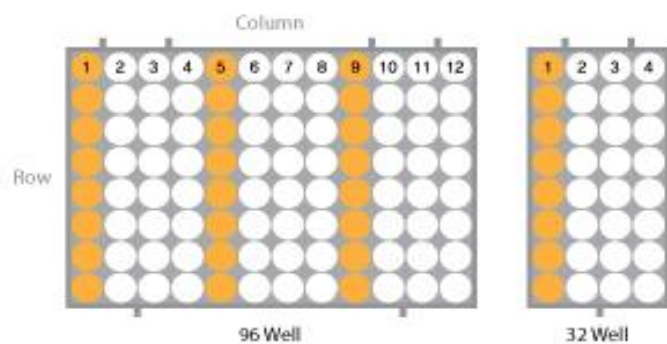


Figure 2: Samples wells (marked in orange) in respect to the Plate type

4. Open the front glass door of the instrument and pull out the plate loader.
5. Insert the rods plastic strip in the strip holder in accordance to the tested samples plate and push it to the end (figure 3).

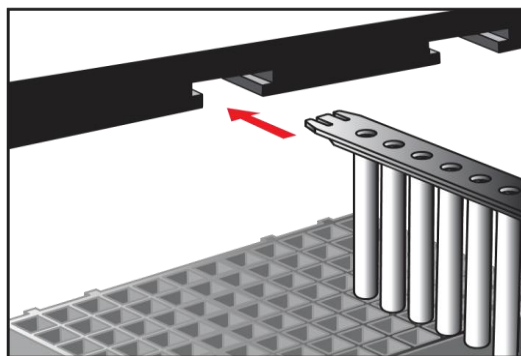


Figure 3: Inserting the rods plastic strip to its position

6. Place the extraction plate on the plate loader. Note that the plastic cavities needs to fit in their position (Figure 4).

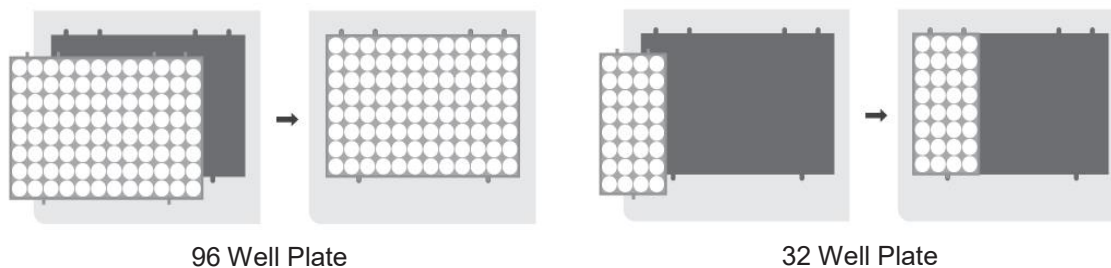


Figure 4: Positioning of the 96 or 32 well plates on the plate loader

7. Push the plate loader in until you hear a clicking sound and pull down the door to close.
8. Select Extraction → Select Round Strip Type → Select Protocol → Select GD → Touch the On / OFF button to select the heating region for the left and right plate loader → Press the Set button to complete the setting → START (Figure 5).

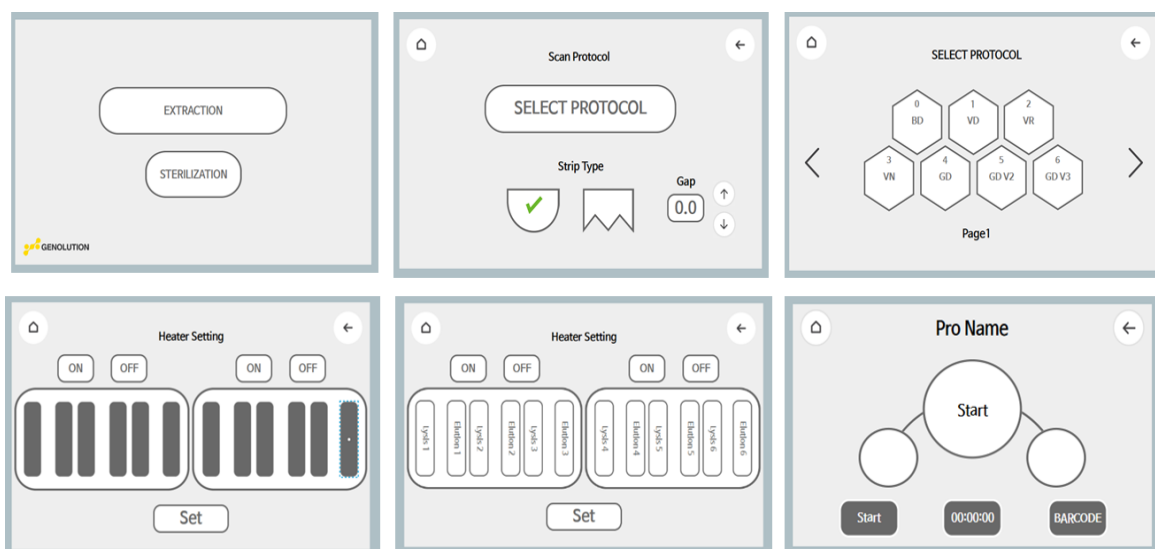


Figure 5: Protocol setting screens

**Attention:** When the Nextractor starts the extraction procedure, it should move the magnetic beads into the wells in the 1<sup>st</sup> column and mix them with the samples.

9. After the extraction is completed, place the Well Plate on Magnetic Bead Holder (provided) and transfer only the eluate to a clean 1.5 mL Centrifuge Tube and proceed to downstream application.

## TROUBLESHOOTING

| Symptoms                                                                                                   | Causes                                                                        | Solution                                                                                                                                                                                                                                                                      |
|------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| When running the extraction protocol the beads are not mixed with the sample in the 1 <sup>st</sup> column | Problem with the instrument calibration position                              | It is recommended to stop the run. Take out the sample plates of the instrument and turn off the instrument power. Then turn on the instrument and start again the run protocol. If the problem reoccur, please refer to the instrument manual for the instrument calibration |
| Considerable loss of elution buffer or no solution found in the well when extraction completed             | Loss of magnetic beads and elution buffer due to excessive gDNA from specimen | It is recommended to reduce the volume of the specimen for extraction and return the run in a new plate                                                                                                                                                                       |
| Excessive magnetic beads found in Elution well                                                             | Including impure material in the specimen                                     | Place the well plate on separated magnetic plate, and use supernatant in the elution well when finished.                                                                                                                                                                      |
| Poor amplification                                                                                         | Carryover Magnetic Beads may interfere with downstream amplification process. | Remove residue Magnetic Beads by using Magnet Bead Holder provided and use supernatant in the Elution Well.                                                                                                                                                                   |
| Cross-contamination                                                                                        | Cross contamination due to user error                                         | Use sterilized laboratory disposables for each sample to prevent sample-to-sample contamination. Avoid splashing when loading the sample into the sample well to minimize contamination of the adjacent wells.                                                                |

For more information, please contact Savyon Diagnostics directly or your Distributor.

[www.savyondx.com](http://www.savyondx.com).

E-mail: [info@savyondx.com](mailto:info@savyondx.com)

## Explanation of Symbols

| Symbols                                                                             | Explanation    | Symbols                                                                             | Explanation                                 |
|-------------------------------------------------------------------------------------|----------------|-------------------------------------------------------------------------------------|---------------------------------------------|
|  | Do not re-use  |  | Please refer to the manual                  |
|  | Lot Number     |  | Expiration Date                             |
|  | CE-IVD marking |  | Manufacturer details                        |
|  | Catalog number |  | Package contains sufficient for < n > tests |
|  | Contents       |  | Authorized European representative          |